

## Go a new way in GC with the BPX90

SGE has introduced the highly polar GC Phase BPX90 consisting of 90 % Cyanopropyl Polysilphenylene-Siloxane, which shows a good stability at a temperature limit of 280°C. Different to „normal“ phases, the BPX90 doesn't separate compounds according to their boiling point but according to their polarity.

This means that double- or triple bonded compounds experience a  $\pi$ - $\pi$  interaction and are stronger separated. Compared to medium polar columns, selectivity is specifically different. It can be expected that this phase separates many compounds differently to what is commonly known which means that the BPX90 can probably solve many individual separation problems.

Interested? Discuss with us if the BPX90 can solve your separation problem! For example at infoExpo 2006.

In addition we would like to invite you to infoExpo in Basel and Lausanne or to Applica, Olten, two upcoming events on analytical technology. Please be aware that the main language of both events is German, except for infoExpo in Lausanne which will be held in French and English.

for further information please refer to

[http://www.infoagil.ch/gc\\_ms/bpx90](http://www.infoagil.ch/gc_ms/bpx90)

<http://www.infoexpo.ch/>


<http://www.sach.ch/applica.html>

# GTG Glastechnik Gräfenroda GmbH

Wiesenweg 35, D-99330 Gräfenroda; Fax +49 36205 77 228; info-Line Tel: +49 36205 77 29; info@gtgvials.de

## Go a new way in GC with the BPX90

The highly polar BPX90 phase separates compounds by their polarity using  $\pi$ - $\pi$  interaction. This allows completely new separations compared to common methods.

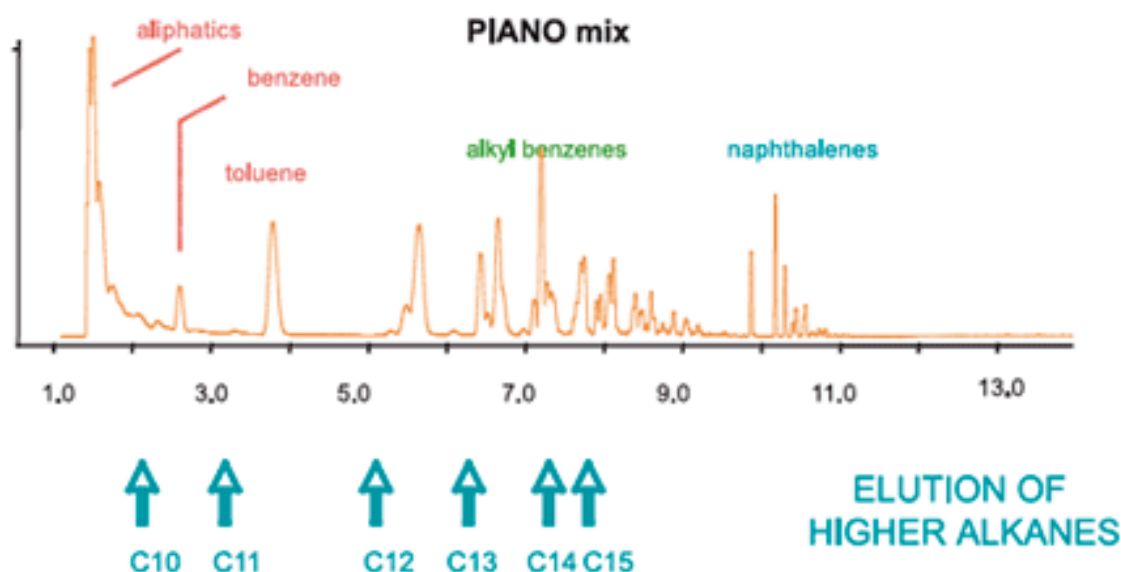

P/N	product nr.	price / pack à
054570	8SBPX90-15/025/025	EUR 360.00 /15 m
054580	8SBPX90-30/025/025	EUR 580.00 /30 m
054590	8SBPX90-60/025/025	EUR 1020.00 /60 m
054573	8SBPX90-15/032/050	EUR 363.00 /15 m
054583	8SBPX90-30/032/050	EUR 649.00 /30 m
054593	8SBPX90-60/032/050	EUR 1120.00 /60 m


all prices are net, excl. VAT  
excl. EUR 10.00 postage and packaging



### APPLICA 2006

to good customers  
we hand out a voucher  
for a free entry -  
value CHF 390.00

conditions apply  
as long as supply lasts

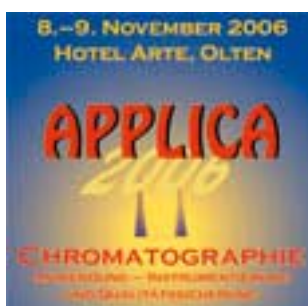

Dienstag  
12. September 2006  
Basel  
UBS Ausbildung- und  
Konferenzzentrum



Mercredi  
13 septembre 2006  
Lausanne  
Hôtel Mövenpick

apply under [www.infomac.ch/applica06](http://www.infomac.ch/applica06)  
more information under [www.sach.ch/applica](http://www.sach.ch/applica)

register under [www.infoexpo.ch/anmeldung](http://www.infoexpo.ch/anmeldung)  
more information under [www.infoexpo.ch](http://www.infoexpo.ch)