

## Robust Macroviai iV2 $\mu$ now also available with screw top

The Macroviai iV2 $\mu$  for very small sample volumes or volatile samples that we introduced about a year ago has proved its worth in the day to day laboratory work. Many critical users now prefer our Macroviai to an insert and vial solution. Therefore we have added to our product line an "HP" screw top version of the iV2 $\mu$ .

Like the  $\mu$ -Vial, the new iV2 $\mu$  Vial has an insert integrated into the body of a 2 ml Vial. In contrast to the  $\mu$ -Vial, however, the insert of the iV2 $\mu$  itself is closed with the cap. This means that the iV2 $\mu$  can be used with volatile samples as no sample will evaporate into the hollow space between vial and insert. The 250  $\mu$ l insert has a very fine and long drawn tip so that the dead volume is less than 5  $\mu$ l.

Our customers especially favour the iV2 $\mu$ -Vial when volatile samples have to be stored in the fridge over an extended period of time or for microreactions where the sample is shaken.



We are happy to send you samples so that you can test the advantages yourself. The iV2 $\mu$  family is now complete. Choose between Crimp, Snap/Crimp or 9 mm Screw Top in clear or amber borosilicate glass. The vials are waiting in stock to be used in your lab.




Soon it's Christmas and we will close our production and office for the holiday. Last working day is Friday 21. December 2007 and first working day in the New Year is 3. January 2008. We thank you - our customers and readers - for your time and your interest in our information and we wish you a happy festive season in the circle of your families and health, happiness and success in the New Year.

# GTG Glastechnik Gräfenroda GmbH

Wiesenweg 35, D-99330 Gräfenroda; Fax +49 36205 77 228; info-Line Tel: +49 36205 77 29; info@gtgvials.de

## iV2 $\mu$ -Vial

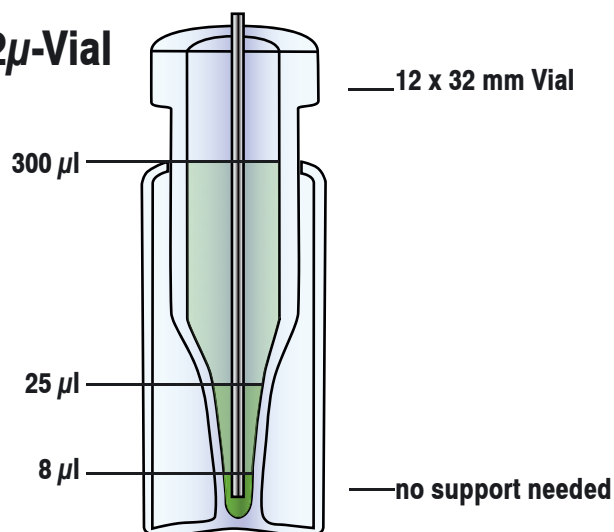

### the advantages at a glance:

#### Insert itself is closed with cap

no sample escapes into the space between vial and insert  
ideal for micro reactions  
suitable for volatile samples and lyophilisation

#### long narrow pointed insert

suited for syringe needles with sidehole  
very low dead volume (< 5 $\mu$ l)

#### established design

stand alone, no support needed  
corresponds with 12 x 32 mm standard  
sturdy design

### Crimp iV2 $\mu$ -Vial



8002-CV-H/iV2 $\mu$   
8002-CV-D/iV2 $\mu$

#### G003-AC-Ru/Te

Aluminium Crimp Cap  
with Rubber/PTFE Septum



### Snap/Crimp iV2 $\mu$ -Vial



8002-SC-H/iV2 $\mu$   
8002-SC-D/iV2 $\mu$

#### 8003-SC-Si/Te

orange Snap Cap  
with Silicone/PTFE Septum



### Agilent compatible iV2 $\mu$ -Vial



8004-HP-H/iV2 $\mu$   
8004-HP-D/iV2 $\mu$

#### 8004-HP-CB-Ru/Te

9 mm blue Screw Cap  
with Rubber/PTFE Septum



Prices in EUR /Price Unit; excl. VAT		Quantity Breaks [pcs.]		
Product No.	Pack of	500 pcs	1000 pcs	10.000 pcs
8002-CV-D/iV2 $\mu$	500 pcs.	57.00 / 100 pcs	54.00 / 100 pcs	51.00 / 100 pcs
8002-CV-H/iV2 $\mu$	500 pcs.	54.00 / 100 pcs	52.00 / 100 pcs	49.00 / 100 pcs
8002-SC-D/iV2 $\mu$	500 pcs.	57.00 / 100 pcs	54.00 / 100 pcs	51.00 / 100 pcs
8002-SC-H/iV2 $\mu$	500 pcs.	54.00 / 100 pcs	52.00 / 100 pcs	49.00 / 100 pcs
8003-SC-Si/Te	100 pcs.	27.00 / 100 pcs	23.10 / 100 pcs	22.20 / 100 pcs
8004-HP-CB-Ru/Te	100 pcs.	22.60 / 100 pcs	19.50 / 100 pcs	18.70 / 100 pcs
8004-HP-D/iV2 $\mu$	500 pcs.	62.00 / 100 pcs	59.00 / 100 pcs	56.00 / 100 pcs
8004-HP-H/iV2 $\mu$	500 pcs.	58.00 / 100 pcs	55.00 / 100 pcs	52.00 / 100 pcs
G003-AC-Ru/Te	500 pcs.	11.60 / 100 pcs	9.90 / 100 pcs	9.50 / 100 pcs

Price excl. VAT, excl. Postage & Packaging